
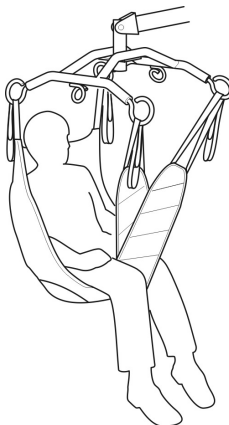


## User Instruction Manual

### Hoyer® Full Back Bariatric Sling

To avoid injury, read user manual prior to use.

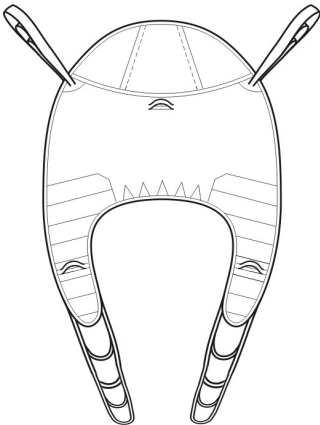

## Introduction

The Full Back Sling is an easy fit, contoured sling designed to fit 85-90% of clients. It is designed for use as the traditional Quickfit Sling but offers greater comfort and negates the need for side suspenders. It has integrated, boned head support and leg padding as standard. The sling has four sizes - medium, large, extra large and XX large.

## Statement of Intended Use

A sling is an item of moving and handling equipment that is used with a mechanical lift in order to facilitate the transfer of a patient. It comprises a specially designed and constructed piece of fabric that is placed under and/or around a patient before being attached to the spreader bar/cradle of a lift to raise, transfer and lower the patient. When selected and used correctly, a sling and lift combination will achieve a safer transfer and reduce the risks associated with manual handling.

It is the responsibility of a competent person to conduct a thorough risk assessment prior to using any sling, to ensure that the sling choice, method of positioning in the sling and procedure for transfer has been correctly determined for the patient. For further guidance, please contact your authorized Hoyer service provider or Joerns Healthcare directly.

**FOR USE WITH PASSIVE LIFTS ONLY.**


## Expected Service Life

The expected service life of a Hoyer sling will vary dependent on use and following the care and washing instructions provided in the user guide. Factors such as wash temperature, detergents, frequency of use and patient weight will impact on the lifetime of your sling. Joerns Healthcare recommend that slings are checked each and every time prior to use to ensure the safety of the patient. Bleached, torn, cut, frayed or broken slings are unsafe and must be discarded and replaced.


## **⚠** WARNING

- **HOYER RECOMMENDS THE USE OF GENUINE HOYER SLINGS.** Hoyer slings and lifters are not designed to be interchangeable with other manufacturer's products. Using other manufacturer's products on Hoyer products is potentially unsafe and could result in serious injury to patient and/or caregiver.
- **For the safety of the patient and carer; before using a sling a full risk assessment must be conducted to ensure that the correct sling choice, method of positioning in the sling and procedure for transfer has been determined for the patient.**
- **CHECK sling and stitching before each use. Using bleached, torn, cut, frayed or broken slings is unsafe and could result in serious injury or death to the patient.**
- **DO NOT** alter slings. Destroy and discard worn slings.
- **NEVER** leave a patient unattended.
- **DO NOT** exceed the rated capacity of the sling or lift.
- **DO NOT** attempt to re-position a patient by pulling on the sling loops.


## How to Fit From a Seated Position

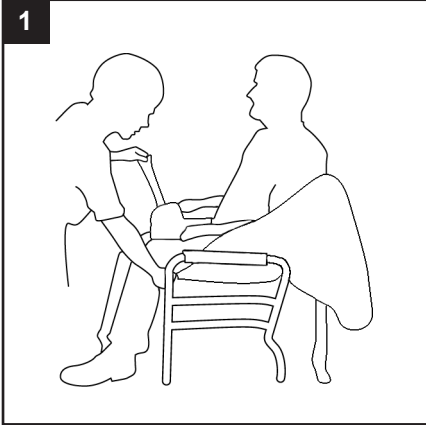

1  
Raise the client's leg and feed the padded leg piece under and up between the legs. Ensure that the sling is not twisted or creased under the thigh. Repeat this procedure for the other leg.

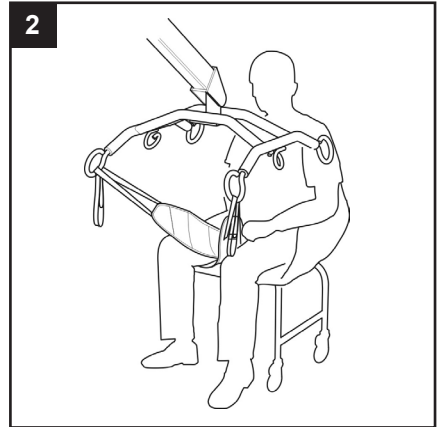

2  
Feed the right leg support through the strap on the left leg support. Cross over leg straps, pass one leg strap through the other and attach to hoist on front hooks.

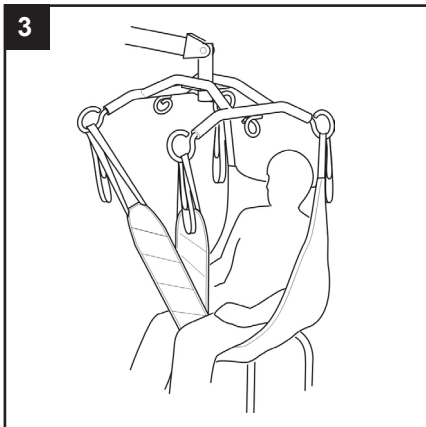


3  
Raise lift to convenient height and attach the shoulder straps as shown. You may then raise the client to the required height. When seating the client the hand grip on the back of the sling may be used to ensure a good seated posture.

### **GOLDEN TIPS**


1. **THE CLIENT'S POSTURE CAN BE ALTERED BY ADJUSTING BOTH SHOULDER STRAPS.**
2. **IF IN DOUBT, ATTACH TO COLOURED LOOPS.**
3. **KEEP SHOULDER STRAPS AT EQUAL LENGTH.**

## How to Fit the Sling from a Lying Position

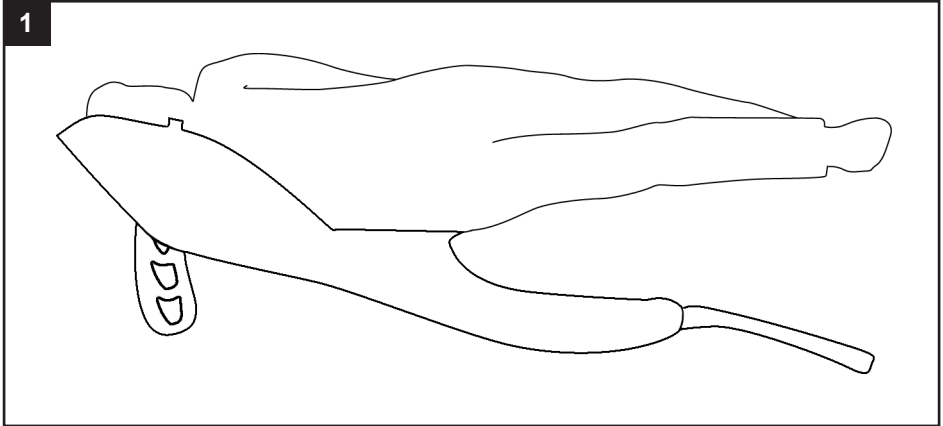

Draw sheet roll the client onto the sling, ensuring that the top of the commode aperture is at the base of the spine.

Bring the leg support straps up and between the client's legs and proceed as from the seated position, attaching to the shortest possible loops. N.B. If you wish to place the client onto a high bed it may be necessary to lower the patient onto an intermediate surface and adjust the strap length.

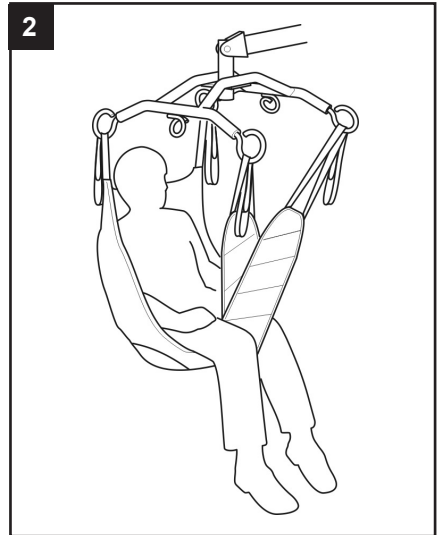

## Technical Specifications

STANDARD	BS EN ISO 10535
----------	-----------------

## Sizing & Safe Working Load

SIZE	FULL BACK BARIATRIC
MEDIUM (YELLOW/WHITE)	850lbs / 386kg
LARGE (GREEN/WHITE)	850lbs / 386kg
EXTRA LARGE (BLUE/WHITE)	850lbs / 386kg
XX LARGE (BLACK/WHITE)	850lbs / 386kg

## Washing Instructions



DO NOT EXCEED 185°F water temperature. Machine wash in warm water.



DO NOT wash with bleach. Bleach will damage the sling's material and make it unsafe for use.



Cool tumble dry, air dry or dry at very low temperature.



DO NOT dry clean.

### **⚠ WARNING**

**Slings can suffer damage during washing and drying and should be checked carefully before each use.**

**PLEASE NOTE:** Additional slings are available to meet individual patient needs. We advise that you always seek the advice of a trained clinician, authorized Hoyer distributor or Joerns Healthcare prior to purchase or use.

THIS PAGE IS INTENTIONALLY BLANK

THIS PAGE IS INTENTIONALLY BLANK

