

Joerns® EasyCare® Series  
leading bed evolution



**Joerns® EasyCare® Series** The new Joerns EasyCare series provides best in class quality and features for a variety of your resident care needs while meeting your facility's budgetary concerns. The EasyCare series adjustable-height low bed has a proven grid-style mattress support platform that provides security for resident positioning, a CenterLine™ steering system for easy maneuvering and a Toroidal frame that is 30 times as strong as standard angle iron steel. The mobility, staff control and assist device/handle options allow facilities to purchase the right bed for the right resident needs, at the right price.

## Features

- 1 Toroidal architecture** engineering of the bed's steel frame is up to 30 times stronger than angle iron
- 2 Therapeutic grid surface** offers improved air flow and ultimate support for the sleep surface
- 3 Four-section mattress support platform** offers full range of comfort and therapeutic positions, which may reduce skin shear
- 4 Attractive bed panel options** available in a wide array of designs and finishes
- 5 Staff control** allows caregivers flexibility to select which features may be adjustable by the resident (option)
- 6 EasyLift Leg Lift** provides quick positioning for improved circulation (EasyCare 5, 7, and 9 models)
- 7 Revolutionary arched legs** deliver a super-strong wide stance for equipment clearance
- 8 Four DC motors** supply amazingly quiet and smooth operation
- 9 7-inch to 30 inch travel range** offers safety, whether resting at an ultra low height or providing care at the higher end (EasyCare 5, 7, and 9 models)
- 10 Assist devices** have superior functionality and perfect grip height
- 11 UltraLock** is standard at the foot end with Centerline steering at the head to maximize security and mobility at any height (EasyCare 7 and 9 only)



EasyCare 7 shown

EasyCare BED SERIES



**EasyCare 9**

- UltraCare XT Base
- EasyCare Deck
- Roll at Any Height
- Optional Staff Control
- UltraLock w/ Centerline steering
- Ergonomic Hand Pendant
- QuickCare



**EasyCare 7**

- UltraCare XT Base
- EasyCare Deck
- Roll at Any Height
- Optional Staff Control
- UltraLock w/ Centerline steering



**EasyCare 5**

- UltraCare XT Base
- EasyCare Deck
- Roll in Low
- Optional Staff Control



**EasyCare 3**


- UltraCare XT Base
- EasyCare Deck
- Roll in Low
- Travel Range of 8¼" to 27"
- Optional Staff Control

ACCESSORIES



**4-Inch Extender**

The 4" extender is designed to extend the mattress support platform an extra 4 inches. This allows the use of the 84" mattress on an 80" bed, or an 80" mattress on a 76" bed.

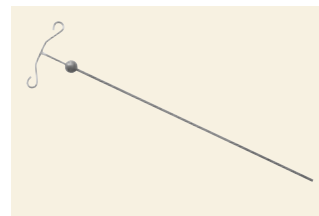

**IV Rod**

The IV rod is adjustable in height to appropriately locate fluids dependent on bed position.



**Floor Mat** The 2" floor mat provides extra security against injury for residents who are prone to falling out of bed.



**Trapeze and Trapeze Adapter**

The trapeze is adjustable to accommodate residents' reach. Additionally, it has three working positions to provide even greater flexibility.

ASSISTS



**Deluxe Assist Handle**

The assist is mounted at an ergonomic height to best assist residents in and out of bed, and has an integrated pendant holder on each side. (EasyCare 9 pendant shown)



**Assist Device**

Features a low profile design, avoiding feeling of restriction or restraint for the resident. Mounted at the perfect ergonomic height to assist residents in and out of bed.



**Assist Handle**

Ergonomically designed to assist residents in and out of bed.

BED SYSTEM OPTIONS



**Joerns PrevaMatt Foam Mattress**

Our Deluxe and Elite mattresses with or without TEMPUR® material feature special mattress construction to keep residents in place when contouring.



**Joerns Air System**

Our complete line of therapeutic air systems include advanced functionality of CairRails and unique design to enhance your wound management program.

# Joerns® EasyCare® Series Specifications

## Frame Construction

- All major structural components, including mattress support platform, are designed with steel tubing and provide the utmost in structural stability

## Additional Information

- All exposed steel surfaces are coated with epoxy powder paint
- Fasteners, screws, rivets and other non-painted parts are corrosion resistant
- All control switches are of momentary contact design
- Head panel is easily lifted off the bed to allow for CPR in emergency situations

## U.S. Patents

- 6473922, 6474869, D445, 614 Other Patents Pending

## Certifications — ETL Listed

- USA:** Conforms to UL 60601-1, 1st ED., 2003
- Canada:** Certified to CAN/CSA C22.2 No. 601.1-M90, CAN/CSA STD C22.2 NO. 60601-2-38-03
- International:** Conforms to IEC 60601-2-38, First ED., 1996, IEC 60601-2-38, Amendment 1, 1999

## Environmental Conditions

### Operating Conditions:

- Ambient Temperature: +10°C to +40°C
- Relative Humidity: 30% to 75% Non-Condensing
- Atmospheric Pressure: 700 hPa to 1060 hPa

### Storage and Shipping Conditions:

- Ambient Temperature: -40°C to +70°C Relative
- Humidity: 10% to 100%
- Atmospheric Pressure: 500 hPa to 1060 hPa

	EasyCare 3		EasyCare 5		EasyCare 7		EasyCare 9	
Bed size	80"	76"	80"	76"	80"	76"	80"	76"
Support surface	35" x 80"	35" x 76"	35" x 80"	35" x 76"	35" x 80"	35" x 76"	35" x 80"	35" x 76"
Overall H length <sup>1</sup>	83½"	80"	83½"	80"	83½"	80"	83½"	80"
Standard overall HF length <sup>1</sup>	85½"	82"	85½"	82"	85½"	82"	85½"	82"
Width <sup>2</sup>	35"		35"		35"		35"	
Low height	8¼"		7"		7¼"		7¼"	
High height	27"		30"		30½"		30½"	
Underbed clearance <sup>3</sup>	19"		23¾"		23¾"		23¾"	
Safe work load <sup>4</sup>	500 lbs.		500 lbs.		500 lbs.		500 lbs.	
Product weight	220 lbs.		220 lbs.		279 lbs.		279 lbs.	
Maximum back angle	65° to horizontal		65° to horizontal		65° to horizontal		65° to horizontal	
Maximum knee angle	35° to horizontal		35° to horizontal		35° to horizontal		35° to horizontal	
Mattress support platform translation - entire travel range	<1"		<1"		<1"		<1"	
Mobility height	8¼"		7"		7¼" - 30"		7¼" - 30"	
Mains voltage	120 VAC		120 VAC		120 VAC		120 VAC	
Frequency	60 Hz		60 Hz		60 Hz		60 Hz	
Input current	3.5 Amps		3.5 Amps		3.5 Amps		4.0 Amps	
Electrical protection	Class I or II, Type B		Class I or II, Type B		Class I or II, Type B		Class I or II, Type B	
Duty cycle	10% (2 min. on/ 18 min. off)		10% (2 min. on/ 18 min. off)		10% (2 min. on/ 18 min. off)		10% (2 min. on/ 18 min. off)	
Circuit protection	Auto reset current sensor		Auto reset current sensor		Auto reset current sensor		Auto reset current sensor	

1 **LENGTH:** "H" is Head panel only; "HF" is Head and Foot panel  
 2 **WIDTH:** With assist device brackets, 36¼"; with assist devices, 41"

3 **CLEARANCE:** Ample clearance between 14" and 30"  
 4 **SAFE WORK LOAD:** Includes bedding, resident, support surface, and all accessories

## Finish

Joerns EasyCare Series comes in a soft tone finish.

We reserve the right to make dimension or specification changes, or deletions consistent with design trends, product improvements, or for reasons beyond our control. Color options have been reproduced as faithfully as printing processes allow.