

BioClinic Therapeutic Support Surfaces

transforming wound management





Imprint[®] Superior support surface designed for use in all care settings

Imprint Features

The ultimate therapeutic foam mattress, the Imprint is made with ultra-soft memory foam that envelops and supports the body while redistributing pressure for continuous patient comfort. This triple-layer, therapeutic foam mattress provides the ideal support surface for prevention and treatment of up to Stage I pressure ulcers.

- Latex free, triple-layer therapeutic foam mattress insures the appropriate foam for each zone of the body, enhancing patient comfort and redistributing pressure
- Incorporates one-inch viscoelastic foam topper to envelop the body and redistribute pressure
- Choice of nylon or stretch cover
- Nylon and stretch covers are treated with an antibacterial, antimicrobial agent to aid with infection control
- Fluid proof and stain resistant cover
- Non-skid bottom
- Meets TB 117 and 16 CFR 1632; with optional fire barrier, meets TB 129, 16 CFR 1633, BFDIX-11
- 350 lb. weight capacity

Flexible Solutions Drive Value

Selecting a support surface that meets the patient's physical capabilities and clinical needs while addressing your financial and risk management issues can be a difficult task. Joerns Healthcare has developed the BioClinic series of five distinct therapeutic foam mattresses to allow you to select the proper surface for both your patient and your budget.

Our BioClinic support surfaces provide:

- Pressure redistribution
- Increased patient comfort and compliance
- A broad selection of surfaces to meet both clinical and financial goals
- Superior craftsmanship for lasting value

Comfort and Safety

Comfort and safety are two primary concerns healthcare providers face when delivering care, and Joerns Healthcare professionals have worked diligently to develop solutions to address these needs when designing our BioClinic series of support surfaces.

Additionally, all BioClinic support surfaces are fully engineered to meet today's stringent regulatory, risk and cost management issues. Joerns Healthcare closely monitors ongoing legislation and takes a proactive stance towards designing our support surfaces to meet current and upcoming regulatory issues.



BodyWrap

Versatile mattress can change along with patient needs

BodyWrap Features

Using a high density, ultra-soft urethane foam that is combined with a firm perimeter, the BodyWrap literally follows its name by surrounding or “wrapping around” the patient for greater comfort and pressure redistribution. By redistributing pressure over a wider surface area, the BodyWrap delivers comfort and is an excellent mattress for patients at risk for pressure ulcers. This true perimeter style mattress provides additional safety, while the removable insert allows for easy replacement to meet the changing needs of the patient.

- Latex free, multi-ply cut and design promotes envelopment and redistributes pressure
- Fully-removable insert to allow caregiver to remove or replace insert to meet the changing needs of the facility or patient
- Firm perimeter-style construction enhances patient safety by reducing the risk of patient falls
- Nylon and stretch covers are treated with an antibacterial, antimicrobial agent to aid with infection control
- Fluid proof and stain resistant cover
- Non-skid bottom
- Meets TB 117 and 16 CFR 1632; with optional fire barrier, meets TB 129, 16 CFR 1633, BFDIX-11
- 350 lb. weight capacity

As a charter member of the Hospital Bed Safety Work Group (HBSW) and an industry leader in entrapment compliance and assistance, we understand the continuous risk and regulatory challenges faced by administrators and staff on a daily basis.

We offer a variety of mattress widths and lengths to accommodate your bed system to help reduce patient entrapment risk.

All Joerns mattresses offer an optional fire barrier in order to comply with Federal 16 CFR 1633 flammability standards.

Quality

Joerns' high quality therapeutic foam mattresses incorporate superior craftsmanship for lasting value. Our full spectrum product range provides solutions for a variety of expenditure levels.

Since the therapeutic benefits of our support surfaces are linked to their ability to bounce back to original shape and size between uses, our foam mattresses are guaranteed not to compress more than 10% during the warranty period (*warranty periods vary by model*).



TurnStyle

Soft, long-lasting, reversible and affordable support surface

TurnStyle Features

Reversible design allows support surface to be flipped or rotated head to toe to extend the life of your support surface, thereby easing your budget. Precision-cut, body contouring waffle cut design provides excellent comfort and eliminates hard spots on the sleeping surface. The result is a comfortable, durable foam mattress with great envelopment and pressure redistribution.

- Latex free, multi-ply therapeutic foam mattress
- Reversible design extends mattress life cycle
- Unique waffle-cut design envelops the body and provides additional patient comfort
- Low shear and low friction cover is treated with an antibacterial, antimicrobial agent to aid with infection control
- Fluid proof and stain resistant cover
- Non-skid bottom
- Meets TB 117 and 16 CFR 1632; with optional fire barrier, meets TB 129, 16 CFR 1633, BFDIX-11
- 350 lb. weight capacity

Wound Management

Advanced Designs, Optimal Outcomes

Since 1976, BioClinic has been a leader in the design of innovative quality mattresses, overlays and positioners for patients in all care settings.

Carefully engineered and tested, Joerns Healthcare's BioClinic products are designed to redistribute pressure and reduce the risk of pressure ulcers. We manufacture our support surfaces with medical grade foam in either a single or multi-ply design utilizing computer-guided, precision-cut production methods. BioClinic support surfaces provide you with quality, comfort, safety, and excellent pressure redistribution to assist you in achieving optimal clinical outcomes.

Customized Solutions

Clinical needs vary. When you purchase a BioClinic mattress, you can customize the size, the cover, and construction style (flat, sidewalls or perimeter). This allows you to provide your standard and bariatric patients with quality, comfort, safety and excellent pressure redistribution.

Our Promise

Support surfaces are an integral component of a complete wound care management program. BioClinic support surfaces provide excellent pressure redistribution and value as part of your wound care program.



BioCore

Great everyday support surface
with advanced heel protection

BioCore Features

An outstanding everyday mattress, the BioCore features our triple-layered construction, incorporating an advanced heel slope zone to enhance comfort and reduce pressure in the delicate heel area.

- Latex free, triple-ply therapeutic support surface
- Precision-cut design promotes envelopment and redistributes pressure
- Incorporated sloped heel section redistributes pressure in delicate heel and lower leg areas
- Available with firm four-inch perimeter style construction to enhance patient safety by reducing the risk of patient falls
- Nylon and stretch covers are treated with an antibacterial, antimicrobial agent to aid with infection control
- Fluid proof and stain resistant cover
- Non-skid bottom
- 350 lb. weight capacity

Construction Style

BioClinic therapeutic foam mattresses are available in three construction styles: flat, sidewall and perimeter. While a flat design is a great choice for most applications, there are some applications that may require a sidewall or perimeter construction.

Sidewall

In the case of the mattress surface being low relative to the bottom of the side rail, the added height of the mattress sidewall will help close the gap. The sidewall also offers added security for the patient when side rails are not used.

Perimeter

This design consists of a high density foam outer rim. Adding the perimeter option may be a good choice for beds that otherwise would fail compliance tests for Entrapment Zone 1. The firmer perimeter resists compression while the center sleeping section provides the patient with therapeutic support.





DermaTech[®] Excellent, everyday value

DermaTech Features

Single-ply, therapeutic high-density foam mattress utilizing precision-cut design, which promotes envelopment and redistributes pressure. Combines quality, comfort and value.

- Latex free, single-ply high-density foam mattress
- Unique stepped contour cuts in leg and heel sections improve pressure reduction in delicate leg and heel areas
- Cover is treated with an antibacterial, antimicrobial agent to aid with infection control
- Soft nylon sewn cover includes a standard 33" foot end zipper (not removable)
- Fluid proof and stain resistant cover
- Non-skid bottom
- Meets TB 117 and 16 CFR 1632; with optional fire barrier, meets TB 129, 16 CFR 1633, BFDIX-11
- 350 lb. weight capacity (6" support surface)
- 250 lb. weight capacity (5" support surface)

Sizing

The right size of mattress is as important as its pressure reducing qualities. By choosing wisely, you can eliminate many of the gaps that plague bed systems.

First, you need an accurate measurement of your bed deck. Then, select a mattress size based on your bed frame. Joerns Healthcare offers Standard, TruFit, and Oversize widths. Standard mattresses are 35" wide. TruFit models are 36" wide and the actual length is ½" longer than your standard bed frames. Oversized

mattresses offer extra surface all the way around the mattress. The Oversize width is a full 36 ¼" and the actual length is 2" longer than your standard bed frames. For example, if you have a 76" bed deck, your mattress will arrive 78". These extra two inches help close the gaps between the bed frame and the end panels.

To help you make the right choice, use the ordering guide on the back page.



Bariatric

Designed for patient comfort
and lasting durability

Bariatric Features

BioClinic bariatric support surfaces are available for the specific needs of the bariatric patient. Our bariatric mattresses are designed with ultra high-density foam, fabricated in layers that will maintain integrity while providing a clinically effective support surface.

- Latex free, designed for patient comfort and lasting durability
- Designed with ultra high-density foam, fabricated in layers to maintain mattress integrity while providing a clinically effective solution
- Low shear and low friction nylon and stretch covers are treated with an antibacterial, antimicrobial agent to aid with infection control
- Meets TB 117 and 16 CFR 1632; with optional fire barrier, meets TB 129, 16 CFR 1633, BFDIX-11
- Available in a variety of sizes to meet the requirements of your bariatric patients
- 800 lb. weight capacity

Covers

Over time, moisture will damage mattress foam. That's why you need covers that will support good wound care and will prevent moisture from leaking through to the foam. BioClinic therapeutic foam mattresses come with a soft, high quality 70-denier nylon cover or an optional transverse stretch cover. Zippered covers are fully replaceable and easily cleaned. All BioClinic nylon mattress covers are low shear and low friction and are treated with an antibacterial, antimicrobial agent to aid with infection control.

Nylon

Our 70-denier covers provide excellent comfort and durability while being fluid proof and stain resistant.

Stretch

The super-soft stretch cover provides greater stretch and conformability as the support surface compresses, making it an excellent choice when pressure management is a challenge and moisture is less of a concern. For additional protection of the support surface, we include a moisture-barrier sleeve around the foam mattress.



BioClinic Series Ordering Information

BioClinic's Foam Mattress Program makes purchasing easy!

Step 1	Choose the model	Factor cost, comfort and pressure redistribution capabilities with the needs of your population
Step 2	Choose the construction style	Consider which style is needed – flat, sidewall or perimeter
Step 3	Choose the width	Determine appropriate width – standard, TruFit or oversize dimension options
Step 4	Choose the length	Consider the beds in your facility and select from 76, 80 and 84 inches
Step 5	Choose the cover	Decide which is best for the application – nylon or transverse stretch fabric

Ordering Information

Model		Construction Style		Width		Bed Frame Length		Cover	
Imprint	IM	Flat	F	Standard	35	76 Inches	76	Nylon Sewn, with End Zipper	NS ¹
BodyWrap	BW	Sidewall	S	TruFit	36	80 Inches	80	Nylon Full Zipper	NZ
TurnStyle	TS	Perimeter	P	Oversize	37	84 Inches	84	Transverse Stretch	SZ ²
BioCore	BC							¹ Available in select styles ² Not available on DermaTech	
DermaTech	DT								

Note: Not all models are available with all construction-styles. See options in chart below.

Sample Order No.


IM S 36 80 NZ BA Optional fire barrier³

model
style
width
length
cover

³All mattresses meet TB 117 and 16 CFR 1632. Meet TB 129, TB 603, 16 CFR 1633, and BFD IX-11 requirements with added fire barrier. Add "BA" to item number to include optional fire barrier. Please note that the BA (fire barrier) option is required to comply with Federal 16 CFR 1633.

Choose Construction Style and Size							Warranties	
Construction Style				Width			Mattress	Warranty Period ⁴
Mattress	Flat	Sidewall	Perimeter	Standard	TruFit	Oversize	Imprint	Ten-years
Imprint	•	•	•	•	•	•	BodyWrap	Ten-years
BodyWrap		•	•	•	•	•	TurnStyle	Ten-years
TurnStyle	•	•	•	•	•	•	BioCore	Five-years
BioCore	•	•		•	•	•	DermaTech	Two-years
DermaTech	•	•		•	•		Bariatric	Two-years

⁴Warranted for defects in workmanship and materials. Please note: Cover warranties are one-year.

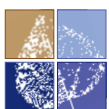

**Joerns
Healthcare**

www.joerns.com • email: info@joerns.com

Post Acute, Acute, HomeCare
5001 Joerns Drive
Stevens Point WI 54481
(P) 800.826.0270
(F) 800.457.8827

VA/Government
7027 Hayvenhurst Avenue
Van Nuys, CA 91406
(P) 800.966.6662
(F) 800.232.9796

Canadian Office
1000 Clarke Road, Suite 6
London, ON Canada N5V 3A9
(P) 866.546.1151
(F) 519.451.8662