



Joerns® P.R.O. Matt Plus Convertible Air Foam Mattress



transforming wound management


Quality. Comfort. Incredibly Advanced.

Our advanced P.R.O. Technology can provide solutions that allow your facility to achieve your clinical and risk management goals.

The P.R.O. Matt[®] Plus is a non-powered mattress replacement system featuring our Pressure Redistribution Optimization (P.R.O.) Technology. With the addition of an optional control unit, the mattress provides powered immersion or alternating pressure therapy, allowing facilities to use one mattress for both pressure ulcer prevention and treatment.

The unique P.R.O. Technology design requires no adjustments or manual inflation. The mattress features a sophisticated internal hosing and valve system that is invisible to the caregiver and patient, but enables the surface to automatically adjust to each patient's weight and body profile.

In the non-powered mode, the P.R.O. Matt Plus is indicated for prevention through Stage II pressure ulcers. With the addition of an optional control unit and based on individual patient assessment, the P.R.O. Matt Plus can be used for the treatment of uncomplicated Stage III or IV pressure ulcers.

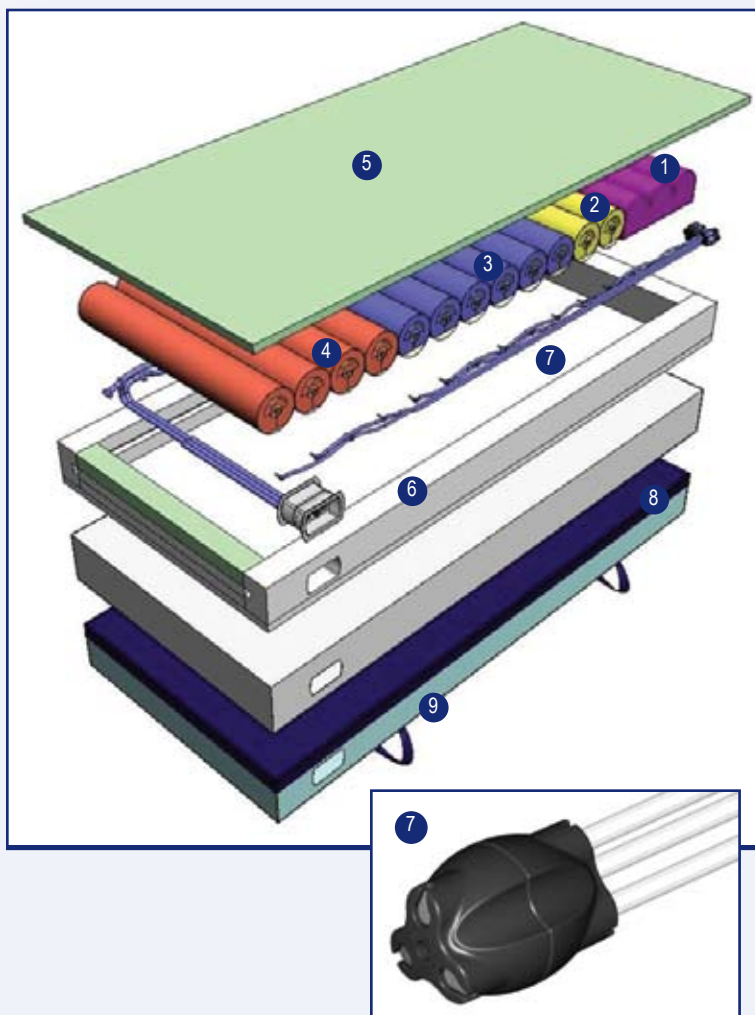

Benefits & Features:

- 1 The head zone is static and comprised of high density foam to optimally support the patient and provide comfort
- 2 The shoulder, torso, and foot zones are comprised of combination air/foam cells. Each zone automatically adjusts to each area of the patient's body to provide optimal pressure redistribution
- 3
- 4
- 5 Viscoelastic foam topper provides maximum pressure redistribution, patient support, and increased comfort
- 6 Surrounding firm perimeter provides stability during patient care and contributes to patient safety by eliminating gaps between frames and siderails
- 7 Custom P.R.O. assembly requires no manual inflation and ensures that each zone remains optimally inflated when the patient is not on the mattress
- 8 Therapy cover is designed to enhance patient comfort and protect the patient's skin:
 - Fluid proof and stain resistant
 - Treated with an antimicrobial to protect it from the growth of mold mildew and bacteria
 - Low shear, low friction nylon cover is standard; stretch cover is optional
- 9 Non-skid, tear resistant bottom prevents mattress sliding on the deck to reduce the risk of falls

Non-Powered P.R.O. Technology

In the non-powered mode, the P.R.O. Matt Plus' internal combination air/foam cells are zoned and designed to reproduce the surface pressure profile of a therapeutic low air loss mattress.

No Maintenance

The P.R.O. Matt Plus automatically re-inflates the air cells in each air zone (shoulder, torso, foot) without the need for caregiver intervention or the use of a pump.

Optional Control Units

By adding an optional control unit, the P.R.O. Matt Plus easily converts into an alternating pressure mattress replacement system. With three control unit options, caregivers can select the unit that best meets their wound management and financial goals.

P.R.O. Matt ES Control Unit Options:

ES and ES2 Control Unit

With the addition of either the ES or ES2 optional control units, the P.R.O. Matt Plus has three modes of operation: *Alternating Pressure*, *Powered* and *Autofirm*.

Alternating Pressure Mode

In the powered mode, the P.R.O. Matt Plus converts to an alternating pressure redistribution system. Three alternating pressure cycle times of 5, 10, and 15 minutes provide customized therapy as needed.

Powered Static Flotation Mode

In this mode, the P.R.O. Matt Plus mattress is not alternating and remains static, allowing for excellent patient immersion and flotation.

Autofirm

To quickly inflate the mattress to maximum pressure, use Autofirm mode. This mode may be desirable to assist with patient transfers and other nursing procedures.

CBC Control Unit

The P.R.O. Matt Plus can also be used in conjunction with the easy-to-use CBC control unit. This pump provides alternating pressure therapy with an adjustable dial. Additional features include static mode and an LED visual alarm alert for low pressure.




Joerns P.R.O. Matt® Plus Convertible Air Mattress

The Joerns P.R.O. Matt Plus dynamic pressure redistribution mattress replacement systems allow healthcare providers to provide optimal interface pressures for low to high-risk patients through controlled air cell inflation. The system easily converts from active alternating pressure redistribution to passive pressure profiling that mimics pressures typically found in a powered therapeutic system.

P.R.O. Matt Plus Ordering Information

Item Number	Description	Cover Option		Zipper	
		Nylon	Stretch	270°	360°
P.R.O. Matt Systems with ES Control Unit					
PMP3680NZBA-ES	36"(W) x 80"(L) x 7"(D)	•			•
PMP3680SZBA-ES	36"(W) x 80"(L) x 7"(D)		•		•
PMP3684NZBA-ES	36"(W) x 84"(L) x 7"(D)	•		•	
PMP3684SZBA-ES	36"(W) x 84"(L) x 7"(D)		•	•	
PMP4280NZBA-ES	42"(W) x 80"(L) x 7"(D)	•		•	
P.R.O. Matt Systems with ES2 Control Unit					
PMP3680NZBA-ES2	36"(W) x 80"(L) x 7"(D)	•			•
PMP3680SZBA-ES2	36"(W) x 80"(L) x 7"(D)		•		•
PMP3684NZBA-ES2	36"(W) x 84"(L) x 7"(D)	•		•	
PMP3684SZBA-ES2	36"(W) x 84"(L) x 7"(D)		•	•	
PMP4280NZBA-ES2	42"(W) x 80"(L) x 7"(D)	•		•	
P.R.O. Matt Systems with CBC Control Unit					
PMP3680NZBA-CBC	36"(W) x 80"(L) x 7"(D)	•			•
PMP3680SZBA-CBC	36"(W) x 80"(L) x 7"(D)		•		•
PMP3684NZBA-CBC	36"(W) x 84"(L) x 7"(D)	•		•	
PMP3684SZBA-CBC	36"(W) x 84"(L) x 7"(D)		•	•	
PMP4280NZBA-CBC	42"(W) x 80"(L) x 7"(D)	•		•	
Mattresses					
PMP3680NSBA	36"(W) x 80"(L) x 7"(D)	•			•
PMP3680SZBA	36"(W) x 80"(L) x 7"(D)		•		•
PMP3684NZBA	36"(W) x 84"(L) x 7"(D)	•		•	
PMP3684SZBA	36"(W) x 84"(L) x 7"(D)		•	•	
PMP4280NZBA	42"(W) x 80"(L) x 7"(D)	•		•	
Control Units					
PMP-ES	ES Control Unit				
PMP-ES2	ES2 Control Unit				
PMP-CBC	CBC Control Unit				

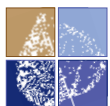

**Joerns
Healthcare**

www.joerns.com • email: info@joerns.com

Post Acute, Acute, HomeCare
5001 Joerns Drive
Stevens Point WI 54481
(P) 800.826.0270
(F) 800.457.8827

VA/Government
7027 Hayvenhurst Avenue
Van Nuys, CA 91406
(P) 800.966.6662
(F) 800.232.9796

Canadian Office
1000 Clarke Road, Suite 6
London, ON Canada N5V 3A9
(P) 866.546.1151
(F) 519.451.8662

© 2009 Joerns Healthcare Inc. • 6220027 RevD • 09-1043