


Joerns
Healthcare

ProCair V • ProCair VR

Therapeutic Mattress Replacement



transforming wound management

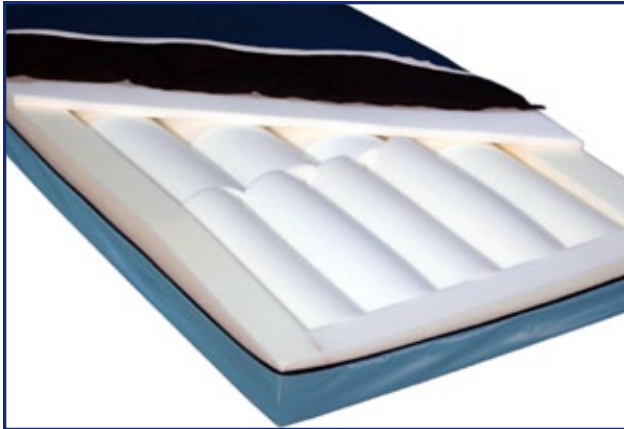

Features and Benefits

Mattress Therapy Features

- Unique longitudinal tubular foam cylinder construction with a high density foam topper for enhanced comfort and pressure redistribution
- The ProCair VR includes a firm perimeter which provides for added stability and safety during patient care and transfer
- Integrated heel slope provides a gentle tapering of the mattress that redistributes pressure along the lower legs while reducing pressure in the delicate heel section
- Ultra-soft Tri-Fill comforter provides added comfort
- Therapy cover with these unique qualities:
 - Fluid proof and stain resistant
 - Low friction and low shear surface
- 70-denier zippered cover includes tie downs and handle straps
- Tear resistant, non-skid bottom prevents mattress sliding on the deck to reduce the risk of falls
- Latex free mattress construction
- Available in multiple widths and lengths to meet individual patient needs
- Psychiatric model available with no zippers, tie-downs or handle straps

Clinical Indications

Recommended for prevention and treatment of Stage I and Stage II pressure ulcers

Technical Specifications

Maximum Weight Capacity*:

400 lbs.

Mattress

Dimensions:

36"(W) x 80"(L) x 8"(D)

36"(W) x 84"(L) x 8"(D)

42"(W) x 80"(L) x 8"(D)

Zipper Cover

360° zipper on 80" mattress

270° zipper on 84" mattress

Warranty

Mattress: Five-year, Cover: One-year

Flammability Standards

Fire barrier meets fire standards TB 117, TB 129, 16 CFR 1633 and BFD IX-11

Density and ILD

Tubular foam cylinders– 1.7 lb. / 20 ILD

Perimeter – 1.7 lb. / 55 ILD

Base Layer – 1.7 lb. / 30 ILD

*Mattress weight capacity only; total weight must not exceed bed frame manufacturer's specified load capacity.

Joerns reserves the right to change specifications without notice.