




Joerns® PrevaMatt® Series
Therapeutic Support Surfaces
transforming wound management



Innovative products creating a safe, comfortable healing environment

The Joerns® PrevaMatt® series of foam mattresses provides residents with comfort and safety while reducing the risk of pressure ulcers. The broad range of Joerns pressure reduction products allows you to choose the right solution to meet clinical, financial and risk management requirements.

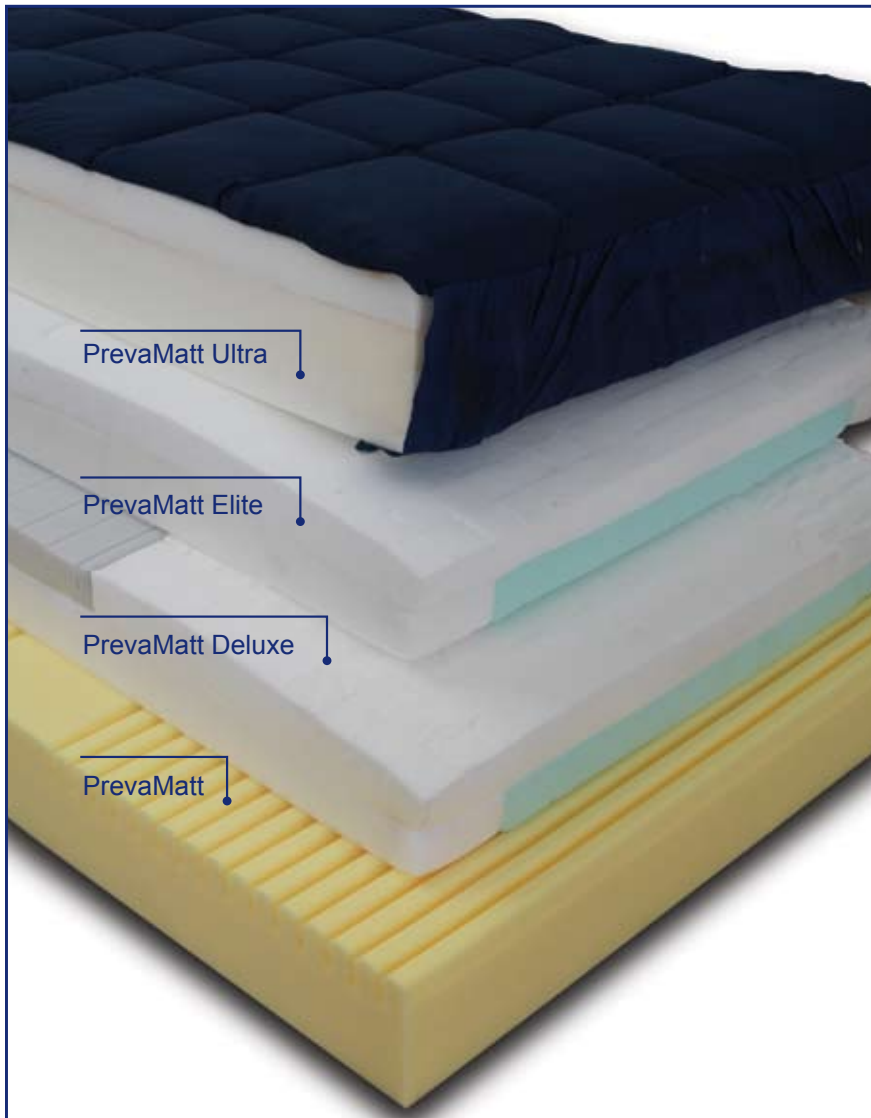

Features:

- Articulation cuts on the Elite and Deluxe models allow the mattress to move in harmony with the bed frame reducing friction and shear on resident skin and increasing product life cycle
- Integrated heel slope provides a gentle tapering of the mattress that redistributes pressure along the lower legs while reducing pressure in the delicate heel section
- 70-denier nylon therapy cover with these qualities:
 - Fluid proof and stain resistant
 - Treated with an antimicrobial to inhibit the growth of MRSA and other bacteria
 - Low friction and low shear surface
 - Non-skid bottom prevents mattress sliding on the deck to reduce the risk of falls, is tear resistant
- Latex free mattress construction
- Fire barrier meets fire standards TB 117, TB 129, 16 CFR 1633 and BFD IX-11

Addressing Resident Safety

Our Joerns surfaces are designed to complete a safer bed system when used in combination with a Joerns bed frame. When used as part of a Joerns bed system, our PrevaMatt mattresses are guaranteed to meet the guidelines set forth in the *Hospital Bed System Dimensional and Assessment Guidance to Reduce Entrapment* developed by the Food and Drug Administration (FDA) and the Hospital Bed Safety Workgroup (HBSW).

To learn more about the FDA entrapment guidelines visit: www.fda.gov/cdrh/beds.

Benefits of our TEMPUR Option


Joerns Healthcare is proud to partner with Tempur-Pedic® to exclusively offer genuine TEMPUR material on our PrevaMatt Elite foam mattresses.

Based on Tempur-Pedic clinical studies, mattresses with made with TEMPUR material are clinically demonstrated to reduce wound incidence rates to less than 1% and to treat up to Stage II pressure ulcers.



PrevaMatt Ultra

- Recommended for prevention and treatment of pressure ulcers up to Stage II
- Unique longitudinal tubular foam cylinder construction with a high density foam topper for enhanced comfort and pressure redistribution
- Ultra-soft Tri-Fill mattress topper provides superior resident comfort and pressure relief
- Firm perimeter provides added stability during resident care or transfer and helps support resident safety
- 400 lb. maximum weight capacity

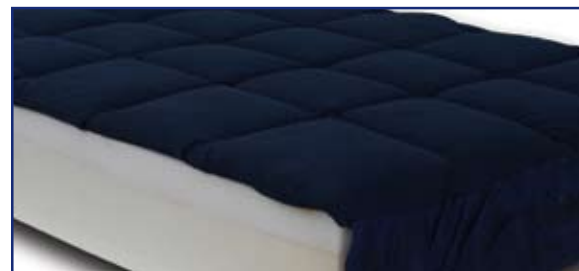

PrevaMatt Elite

- Recommended for prevention and treatment of pressure ulcers up to Stage II
- Two-ply construction Incorporates cross-cut two-inch of viscoelastic foam topper to envelope the body and redistribute pressure
- Optional TEMPUR material available for entire mattress topper
- Firm perimeter and Tru-fit sizing help meet FDA/HBSW guidelines
- Strategically located mattress base cuts promote easier flexing during bed frame articulation - reduces wear and tear to help extend mattress life
- 500 lb. maximum weight capacity



PrevaMatt Deluxe

- Recommended for prevention and treatment of pressure ulcers up to Stage I
- Multiple-ply therapeutic high-density foam mattress with soothing, viscoelastic foam in heel slope for delicate heel section.
- Firm perimeter and Tru-fit sizing help meet FDA/HBSW guidelines
- Strategically located mattress base cuts promote easier flexing during bed frame articulation - reduces wear and tear to help extend mattress life
- 500 lb. maximum weight capacity



PrevaMatt

- Recommended for prevention and residents at risk of developing pressure ulcers
- Single-ply, therapeutic high-density foam mattress which encourages envelopment and redistributes pressure
- 300 lb. maximum weight capacity



Added Security with Sidewalls

In the case of the mattress surface being low relative to the bottom of the side rail, the added height of the mattress sidewall will help close the gap and help to reduce the risk of entrapment.

The sidewall can also aid in fall prevention by providing added security for the patient when side rails are not used.



Firm Perimeter Provides Safety

This design consists of a high density foam outer rim. Adding the firm perimeter option may be a good choice for beds that otherwise would fail compliance tests for the FDA Guidelines Entrapment Zone 1.

This firm perimeter can help reduce the risk of falls by providing stability as the resident sits on the edge of the bed or during ingress and egress.



Joerns Bed Systems A Single Solution for your Resident

Joerns Healthcare designs our healthcare beds, support surfaces and accessories to be used as a complete bed system - the proper combination of frame, support surface, head/foot board and assist devices are needed to minimize the risk of entrapment and comply with the *2006 FDA Hospital Bed Dimensional and Assessment Guidance to Reduce Entrapment*.

We are very confident in our bed systems. For this reason, when you purchase a Joerns bed system, we guarantee compliance with these entrapment guidelines. Call 800.826.0270 for details.



Benefits of a Joerns Bed System

- Broad product selection gives you the flexibility to customize the mattress to create the solution that will best address your resident's clinical and comfort needs
- Optional mattress TruFit sizing reduces the risk of life-threatening entrapment issues
- Hi-Low bed feature helps prevent resident falls
- Manufactured in the USA to meet FDA, OSHA, ISO, IEC, UL, and additional self-imposed standards that all promote the quality and functionality of our products leading to overwhelming customer satisfaction
- Single service solution provides value added-services including clinical and technical solution-based programs to help ensure your resident is safe and your facility is protected

Introducing the Evolutionary UltraCare® XT Bed

We already make the best bed in the business. Naturally, we wanted to improve it and increased functionality was our goal. The new UltraCare XT is safer, more comfortable, faster and easier to use, with advanced configurability features no other bed offers. We've added to the great features you've trusted in the UltraCare bed.



AdvanceCare – Features gravity-assist for easier and safer repositioning. Utilize Comfort Chair functionality for comfort, relaxation, improved circulation, and more ergonomic care by bringing the resident closer to the caregiver.

UltraWide – Extend the resident's sleep width to 39" (twin size) or 42" without tools! Residents love the extra space for comfort. Your staff will love it because it helps to deliver safer transfers and helps to prevent falls throughout the night. You can even accommodate larger residents, reducing the need for bariatric equipment!

UltraLock – Designed for mobility or security at any height, while allowing for a 7" low height. UltraLock provides quick, easy bed movement.

Enhanced controls – Enhanced staff control design feature IPX4 fluid ingress protection rating, angled view and optional backlighting for easy operation day or night.

Residential-Style wood trim kit – Attract new residents with a room décor that rivals that of home. No tools addition and removal of wood trim to your UltraCare XT frame.

UltraCare XT: It's One Bed. A Total Solution. Call Today 800.826.0270.

