

Joerns
Healthcare

ReplaMatt Basic

Therapeutic Mattress Replacement



transforming wound management



Features and Benefits

Mattress Therapy Features

- Firm, yet highly resilient foam, gently cradles head and torso for maximum patient comfort
- Integrated viscoelastic foam heel slope provides a gentle tapering of the mattress that redistributes pressure along the lower legs while reducing pressure in the delicate heel section
- Strategically located articulation cuts in mattress base allow the mattress to move in harmony with the bed frame reducing friction and shear on patient skin and increasing product life cycle
- Firm perimeter and Tru-fit sizing help meet FDA/HBSW guidelines
- 70-denier therapy cover with these unique qualities:
 - Fluid proof and stain resistant
 - Infection control treatment
 - Low friction and low shear surface
- Non-skid bottom prevents mattress sliding on the deck to reduce the risk of falls, is tear resistant
- Latex free mattress construction

Clinical Indications

Pressure Relief: Recommended for prevention and treatment of pressure ulcers up to Stage I.

Technical Specifications

Maximum Weight Capacity*: 500 lbs

Mattress

Dimensions: 35"(W) x 80"(L) x 7"(D)

36"(W) x 80"(L) x 7"(D)

36"(W) x 84"(L) x 7"(D)

Warranty

Mattress: Five-years

Cover: One-year

Flammability Standards

Fire barrier meets fire standards TB 117, TB 129,
16 CFR 1633 and BFD IX-11

ILD and Density


Top Layer – 1.7 lb. / 15 ILD

Top Layer Heel – 3.0 lb. / 15 ILD

Perimeter – 1.7 lb. / 55 ILD

Base Layer – 1.7 lb. / 30 ILD

**Mattress weight capacity only; total weight must not exceed bed frame manufacturer's specified load capacity.
Joerns reserves the right to change specifications without notice.*

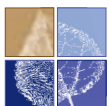

**Joerns
Healthcare**

www.joerns.com • email: info@joerns.com

Post Acute, Acute, HomeCare
5001 Joerns Drive
Stevens Point WI 54481
(P) 800.826.0270
(F) 800.457.8827

VA/Government
7027 Hayvenhurst Avenue
Van Nuys, CA 91406
(P) 800.966.6662
(F) 800.232.9796

Canadian Office
1000 Clarke Road, Suite 6
London, ON Canada N5V 3A9
(P) 866.546.1151
(F) 519.451.8662

© 2009 Joerns Healthcare Inc. • 6220192 RevA • 09-0458