


Expandable Bariatric Bed Frame

As the brawnier version of our RC 750, the RC 1000 is an economical 1000 lb capacity bed frame with five function electric control. The simple to use pull pin system allows caregivers to easily expand the bed frame's width from 39" to 54" to accommodate various patient girth and emergency transport through 40" door frames. Its split frame design allows for easy delivery into tight areas. The RC 1000's reasonable price makes it an excellent bariatric bed solution for the alternate care market.



Multi-functional Bariatric Capacity in an Affordable Package

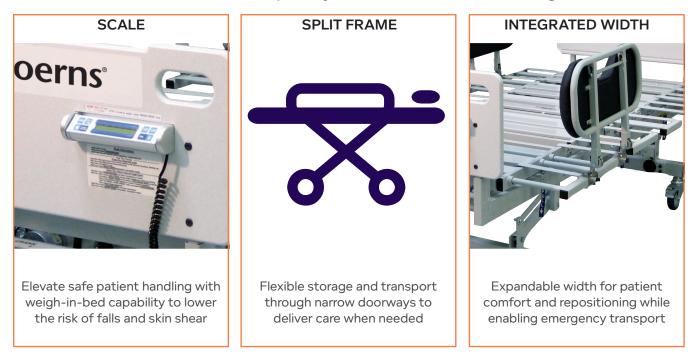



Standard Features

- Split frame design for easy delivery and storage
- 5 functions including Trendelenburg and cardiac chair
- Integrated scale for one-button weighing in any position
- Expandable bed deck permits transport through doorways as narrow as 40"
- Accommodates mattresses up to 54" wide
- 3 position, wing type, padded side rails
- IV pole and Foley bag holder
- 1000 pound weight capacity
- Optional trapeze

Frame Specifications

	- ···	
Bed Deck Length:	84"	
	(96" with optional	extender)
Overall Frame Size:	39.9" to 56.5" W to 91" L	
Mattress Size:	Minimum	37"
	Maximum	54"
Bed Deck Height:	Minimum	17"
	Maximum	30"
Head/Foot Panel Width:		38"
Safe Working Load:		1000 lbs
Power Requirements: 115 VAC, 60 Hz, 10.5 AMP		
Warranty:	Structural	5 years
	Electrical System	1 year

.

RC 1000 Suite Package

Includes:

- RC1000
- Pressure redistribution mattress with Safety Perimeter
- Bed Alarm for patient exit notifications
- Fall pad for added piece of mind

UPGRADE THE SUITE with a lift or CareChair for safe patient handling





