

## Versatile and Portable Stand Aid with Compact Footprint to Allow Navigation in Confined Spaces

The Journey is an ultra-compact, folding stand aid with best-in-class functionality. Both lightweight and nimble, it can be folded for easy storage and transportation. The unique adjustable cradle supports a wider range of patient heights and sizes, allowing your facility to have one lift that meets the need of many.



### UNIQUE FOLDING DESIGN



Compact folding design enables transportation outside the confines of the home and allows the Journey to stand unassisted when stored.

### ADJUSTABLE CRADLE



Unique adjustable cradle supports a range of patient heights and sizes and lends itself to transfers from lower surfaces.

### ADJUSTABLE KNEE PAD



Contoured, height adjustable knee pad offers high level comfort and accommodates a large range of patient sizes.


### REMOVABLE FOOT REST



The foot rest transforms the lift into a valuable restorative care aid, allowing the patient to fully stand and walk under close supervision.



**SAFE WORKING LOAD: 340 LB / 155 KG**

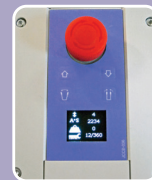

### HOYER SMART MONITOR

- ◊ Lift Cycle Counter
- ◊ Overload Protection
- ◊ Service Reminder
- ◊ Emergency

### SPECIFICATION

Maximum Overall Length	39.4"	1000 mm
Minimum Overall Length	36.8"	935 mm
Maximum Overall Height	63"	1600 mm
Minimum Overall Height	44.5"	1130 mm
Legs Open - External Width	38.2"	970 mm
Legs Open - Internal Width	34.5"	875 mm
Legs Closed - External Width	27"	685 mm
Legs Closed - Internal Width	22.2"	590 mm

### WEIGHTS

Mast, Base and Boom Assembly	82.4 lb	37.4 kg
Power Pack	6.2 lb	2.8 kg
Total	88.6 lb	40.2 kg
Base Assembly	30.2 lb	13.7 kg
Mast & Boom Assembly (no battery)	42.3 lb	19.2 kg

Specifications may vary due to manufacturing tolerances.