

PERSONALIZED COMFORT



Personalized AutoTransfer Height promotes resident dignity and facilitates safe mobility, including bariatric individuals.



Unique 30°/45° comfort pauses reminders to aid sleep, circulation and respiratory therapy.



Adaptable length extensions and width expansions serve more residents.


VERSATILITY



Hand pendant with integrated hook meets the challenges of institutional care with tactile, back-lit buttons for darkened rooms.



Sleep Mode feature saves 80% more power when bed is not in use.

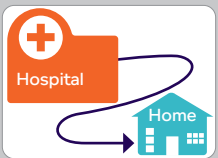

Swing bed flexibility tailors care for everyday residents plus dignified bariatric care transitioning hospital to home.

PROTECTION



Fast, care-assisted repositioning to aid higher acuity admissions and encourage continuous mobility.



Patented arched leg⁺ offers industry-leading low bed height that mitigates fall injury.



LED SafeLight^{*} provides floor visibility in three soft intensities for safe bed egress at night.



"ULTRA" MEANS UNRIVALLED, BEST-IN-CLASS:
 -WEIGHT LOAD (SWL 650 LBS)
 -TRAVEL RANGE
 -FRAME STRENGTH
 -STABILITY



SAFE WORKING LOAD: 650 LB

SPECIFICATION SNAPSHOT	
Frame Mobility	Roll-in-Low (RIL) Roll-any-Height (RAH)
Height from Floor	7" - 30"
Dimensions (W x L)	35" x 80"
Overall Bed Length	85.5" L
Pendant Functions	Back with Comfort Pause, knee, hi-low, Comfort Chair, AutoContour, personalized AutoTransfer, lockout with indicators
Staff Control Functions*	Pendant functions and Gravity Assist, Trendelenburg, reverse Trendelenburg, personalized AutoTransfer Height, function lockout, SafeLight Control
Sleep Deck Style	Toroidal steel therapeutic grid
Actuators	4 actuator system
Casters	Heavy-duty 3" diameter
Horizontal Movement	< 1"
Battery Backup	Built-in*
WARRANTY INFORMATION (see user manual)	
Electrical Components	5-years
Bed Frame	15-years
Welds	Lifetime

Joerns[®] Connexio[®] smart network enables improved efficiency and delivery of care by empowering operators and clinicians with real time visibility of assets, staff and residents. Ready to connect your EasyCare[®] and UltraCare[®] beds?



Connexio[®]
Smart Network
 Care Management,
 Simplified

FRAME	
Frame Mobility	Roll-in-Low (RIL), Roll-any-Height (RAH)
Mobility Lock Mechanism	RIL: Large, High-grip Rubber Feet + Locking Casters, RAH: UltraLock® Three-position Mobility Lock
Overall Bed Length	85.5", 89.5", 93.5"
Safe Working Load	650 lb
Product Weight	RIL: 220 lb, RAH: 247 lb/274 lb
SLEEP DECK	
Height from Floor	7" - 30"
Dimensions (W x L)	35" x 80"
Width Expansion Options	Fixed Width Accessory: 42"
Length Extension Options	4", 8"
Angle Range	Back: 0 - 65° (79° max. in ComfortChair), Knee: 0 - 35°
ELECTRICAL	
Power Saving Standby Mode	0.5W at idle
Circuit Protection	Auto Reset Current Sensor
Backup Battery Option*	On-board Charging
ADVANCECARE POSITIONING	
Standard with Backlit Control	Hi-Low, Back, Knee, One-touch Comfort Chair, AutoContour, One-touch Programmable AutoTransfer, Programmable Low Height, Pendant Lockout Access with Indicators
Back Angle and Low Height Pauses	Standard Comfort Pauses (30° and 45°), Low Height Obstruction Reminder
Included with Staff Control	AdvanceCare Staff Control with Function Lockout, Programmable AutoTransfer Height, Programmable Gravity Assist (4°) and Trendelenburg (12°)/Reverse Trendelenburg(12°), SafeLight Under-bed Light Control*, 3 level programmable
COMPLIANCE	
HBSW Compliant	Bed System
Regulation Compliance	Conforms to IEC 60601-1, 60601-1-2, and 60601-2-52

Specifications may vary due to manufacturing tolerances.

* Sold Separately

† US Patent #6473922, #D437143

ASSIST DEVICES



F028 Deluxe Assist Handle



F026 Pivoting Assist Handle



F17 Assist Device



Trapeze Patient Helper

ACCESSORIES



Multiple Panel Style and Finish Options



AdvanceCare Staff Control




N641 LED SafeLight



N111 Battery Backup



42" Width Expansion

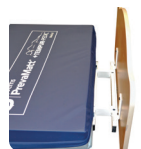

4" / 8" Length Extension